



Neogen Corporation
944 Nandino Blvd., Lexington KY 40511 USA
800/477-8201 USA/Canada | 859/254-1221
Fax: 859/255-5532 | E-mail: inform@neogen.com | Web: neogen.com

COTININE

ELISA KIT INSTRUCTIONS PRODUCT #182219, 182215 & 182213
FORENSIC USE ONLY

INTENDED USE: For the determination of trace quantities of Cotinine and/or other metabolites in human serum and other forensic samples.

DESCRIPTION

Neogen Corporation's Cotinine ELISA (Enzyme-Linked ImmunoSorbent Assay) test kit is a qualitative kit designed for use as a screening device for the detection of drugs and/or their metabolites. The kit was designed for screening purposes and is intended for forensic use only. It is recommended that all suspect samples be confirmed by a quantitative method such as gas chromatography/mass spectrometry (GC/MS).

ASSAY PRINCIPLES

Neogen Corporation's test kit operates on the basis of competition between the drug or its metabolite in the sample and the drug-enzyme conjugate for a limited number of antibody binding sites. First, the sample or control is added to the microplate. Next, the drug-enzyme conjugate is added and the mixture is incubated at room temperature. During this incubation, the drug in the sample or the drug-enzyme conjugate binds to antibody immobilized in the microplate wells. After incubation, the plate is washed to remove any unbound sample or drug-enzyme conjugate. The presence of bound drug-enzyme conjugate is recognized by the addition of K-Blue® Substrate (TMB). After a 30 minute substrate incubation, the reaction is halted with the addition of an acid stop. The test can be read visually or with a microplate reader equipped with a 450 nm filter. The extent of color development is inversely proportional to the amount of drug in the sample or control. In other words, the absence of the drug in the sample will result in a dark yellow color, whereas the presence of the drug will result in light yellow to no color development.

STORAGE AND STABILITY

This kit can be used until the expiration date on the label when stored refrigerated at 2-8°C.

MATERIALS PROVIDED – SINGLE KIT (96 WELL)

1. **Wash Buffer Concentrate (10X):** 20 mL. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
2. **K-Blue Substrate:** 20 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H₂O₂) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
3. **Drug-Enzyme Conjugate:** 200 µL. Drug-horseradish peroxidase conjugate. Dilute 100X before use.
4. **Drug-Enzyme Conjugate Diluent:** 15 mL (ready-to-use). Provided for dilution of drug-enzyme conjugate.
Note: Diluent is specific for the Cotinine drug-enzyme conjugate.
5. **Antibody Coated Plate:** A 96-well Costar plate, in strips of 8 break-away wells, coated with anti-drug antiserum. The plate is ready for use as is. Do not wash.
6. **Acid Stop Solution:** 18 mL (ready-to-use). 1N H₂SO₄ used to stop the enzyme reaction.
7. **QC Positive Control:** 500 µL provided (synthetic matrix), 25 ng/mL. Do not dilute.
8. **QC Negative Control:** 500 µL provided (synthetic matrix), 0 ng/mL. Do not dilute.

MATERIALS PROVIDED – BULK KIT (480 WELL)

1. **Wash Buffer Concentrate (10X):** 100 mL. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
2. **K-Blue Substrate:** 100 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H_2O_2) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
3. **Drug-Enzyme Conjugate:** 1 mL. Drug-horseradish peroxidase conjugate. Dilute 100X before use.
4. **Drug-Enzyme Conjugate Diluent:** 85 mL (ready-to-use). Provided for dilution of drug-enzyme conjugate.
Note: Diluent is specific for the Cotinine drug-enzyme conjugate.
5. **Antibody Coated Plate:** 5 X 96 well Costar plates, in strips of 8 break-away wells, coated with anti-drug antiserum. The plates are ready for use as is. Do not wash.
6. **Acid Stop Solution:** 90 mL (ready-to-use). 1N H_2SO_4 used to stop the enzyme reaction.
7. **QC Positive Control:** 2 mL provided (synthetic matrix), 25 ng/mL. Do not dilute.
8. **QC Negative Control:** 2 mL provided (synthetic matrix), 0 ng/mL. Do not dilute.

MATERIALS PROVIDED – BULK KIT (4800 WELL)

1. **Wash Buffer Concentrate (10X):** 1 L. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
2. **K-Blue Substrate:** 2 X 500 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H_2O_2) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
3. **Drug-Enzyme Conjugate:** 9.5 mL. Drug-horseradish peroxidase conjugate. Dilute 100X before use.
4. **Drug-Enzyme Conjugate Diluent:** 850 mL (ready-to-use). Provided for dilution of drug-enzyme conjugate.
Note: Diluent is specific for the Cotinine drug-enzyme conjugate.
5. **Antibody Coated Plate:** 50 X 96 well Costar plates, in strips of 8 break-away wells, coated with anti-drug antiserum. The plates are ready for use as is. Do not wash.
6. **Acid Stop Solution:** 2 X 500 mL (ready-to-use). 1N H_2SO_4 used to stop the enzyme reaction.

MATERIALS NEEDED BUT NOT PROVIDED

1. Deionized water.
2. Precision pipettes that range from 10 μ L - 1000 μ L and disposable tips.
3. Graduated cylinder to dilute and mix wash buffer.
4. Plate cover or plastic film to cover plate during incubation.
5. Clean glassware (i.e. test tubes) to dilute samples.
6. Microplate reader with 450 nm filter.
7. Cutoff calibrator.
8. Dilution Buffer (1X PBS) for urine samples.
9. Oral Fluid Buffer for oral fluid samples.
10. EIA Buffer for breast milk samples.

OPTIONAL MATERIALS

1. Microplate shaker

PRECAUTIONS AND NOTES

1. **DO NOT** use kits or components beyond expiration date.
2. **DO NOT** mix conjugates and plates from different kit lots.
3. **DO NOT** pipette reagents by mouth.
4. Pour K-Blue Substrate out of the bottle into a clean reservoir. To prevent contamination of the substrate, **DO NOT** pipette out of the bottle.
5. All specimens should be considered potentially infectious. Exercise proper handling precautions.
6. Keep plate covered except when adding reagents, washing or reading.
7. Kit components should be refrigerated at all times when not in use.
8. Use aseptic technique when opening and removing reagents from vials and bottles.
9. **DO NOT** smoke, eat or drink in areas where specimens or reagents are being handled.
10. **DO NOT** substitute DI water for the wash step of this protocol. Use only Neogen's wash buffer.

11. Sodium Azide concentrations at 0.01% or less should not interfere with the assay provided that recommended dilutions are followed.
12. **DO NOT** reuse wells, they are for one use only.

PROCEDURAL NOTES

1. Desiccant bag must remain in ziplock pouch with unused strips. Keep ziplock pouch sealed when not in use to maintain a dry environment.
2. Use clean pipette tips for the buffer, drug-enzyme conjugate, controls and samples.
3. Before pipetting a reagent, rinse the pipette tip three times with that reagent.
4. When pipetting into the wells, **DO NOT** allow the pipette tip to touch the inside of the well or any of the reagent already inside the well. This may result in cross contamination.
5. Controls and samples should be assayed in duplicate.
6. Before substrate addition, wipe the outside bottom of the wells with a lint-free wiper to remove dust and fingerprints.
7. Gently mix specimens and reagents before use. Avoid vigorous agitation.

SAMPLE TREATMENT

Recommended sample preparations are listed below. These may change based on your laboratory's determination.

- a. **Serum:** Can be assayed directly, undiluted.
- b. **Human Breast Milk:** Dilute 1:5 in EIA Buffer (available for purchase from Neogen Corporation as item number 301276 (500mL) or 301277 (1L)).
- c. **Oral Fluid:** Dilute 1:4 in Oral Fluid Buffer (available for purchase from Neogen Corporation as item number 120176 (500mL) or 120177 (1L)).
- d. **Urine:** Dilute human urine samples 1:50 in Dilution Buffer (1X PBS); available for purchase from Neogen Corporation as item number 27029).
- a) **Other Forensic sample types:** Please contact your Neogen Representative for assistance.

TEST PROCEDURES


The following test procedures can be run manually or on an automated instrument. Please contact your Neogen representative for assistance with protocols for automated instruments.

1. Prior to performing the test, allow 60 minutes for all kit components and samples to reach room temperature, i.e., 20-23°C (68-74°F). Mix well by gentle inversion prior to use.
2. Determine the number of wells to be used.
3. Prepare the enzyme conjugate solution by diluting the enzyme conjugate stock 1 to 100 in the Cotinine Conjugate Diluent provided. Mix the solution by inversion. Do not vortex. Example: for two eight well strips, add 25 µL of the enzyme conjugate stock to 2475 µL of Drug-Enzyme Conjugate Diluent. Store unused conjugate at 2-8°C. For convenience, fully diluted enzyme can be prepared up to 5 days prior to use.
4. Add 20 µL of sample, Neogen controls or laboratory calibrators to the appropriate wells in duplicate. **DO NOT** dilute Neogen's positive or negative controls
5. Add 100 µL of diluted drug-enzyme conjugate to each well. Use 8-channel pipetter or 12-channel pipetter for rapid addition.
6. For manual runs, mix by gently shaking plate. A microplate shaker may be used.
7. Cover plate with plastic film or plate cover and incubate at room temperature for 30 minutes.
8. During the conjugate incubation, dilute concentrated wash buffer 10 fold with deionized water (i.e. 20 mL of concentrated wash buffer plus 180 mL of deionized water). Mix thoroughly. Diluted wash buffer is stable for 5 days at room temperature or 7 days at 2-8°C.
9. Once the incubation is complete, dump or aspirate the liquid from the wells. Tap the plate on a clean lint-free towel to remove any remaining liquid in the wells.
10. Wash each well with 300 µL of diluted wash buffer. Manual Wash: For manual wash procedures repeat for a total of 3 washings, invert and tap dry the plate following each step. After completing the last wash step wipe the bottom of the wells with a lint-free towel to remove any liquid on the outside of the wells. Automated Wash: If an automated plate washer is used wash the plate for a total of 5 washings with 300 µL of diluted wash buffer. It is important for the automated washer to conduct a final aspirate cycle to eliminate residual amounts of wash buffer. Residual amounts of buffer in the wells will affect assay performance. Note: DI water should never be used for the plate wash.
11. Add 150 µL of the K-Blue Substrate to each well. For manual runs, use a multi-channel pipetter for best results.
12. Incubate at room temperature for:

Serum	15 Minutes
Human Breast Milk	15 Minutes
Oral Fluid	10 Minutes
Urine	8 Minutes

13. Add 150 μ L of the Acid Stop (1N H₂SO₄) to each well to stop enzyme reaction. Mix gently before measuring the absorbance. For automated systems a 10 second shake is sufficient. Measure the absorbance at a wavelength of 450 nm. Wells should be read within 2 hours of stopping the reaction.

STANDARD CURVE IN SYNTHETIC MATRIX

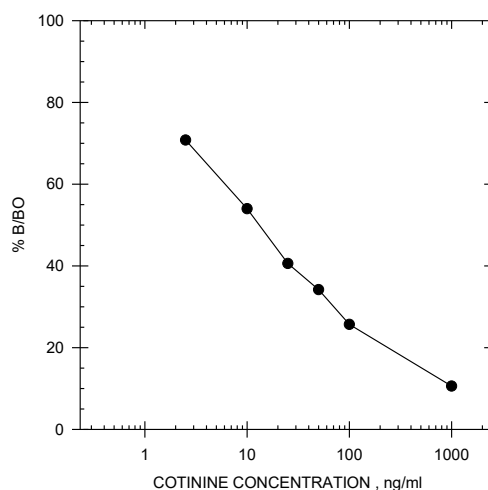

SENSITIVITY

Compound	I-50 in Synthetic Matrix	I-50 in Human Breast Milk	I-50 in Human Urine	I-50 in Oral Fluid
(-)Cotinine	39.4 ng/mL	23.0 ng/mL	10.2 ng/mL	12.5 ng/mL

The term I-50 is used to define the sensitivity of the test. This number is derived from a standard curve generated with the drug in synthetic matrix. The drug concentration that shows 50% less color activity than the zero standard is considered to be the I-50.

SPECIFICITY

Compound	Compound Concentration (ng/mL)	(-)Cotinine Equivalents (ng/mL)	% Cross Reactivity
(-) Cotinine	39.4	39.4	100
(+/-) Trans-3'-hydroxycotinine	633	39.4	6.5
(R/S) Norcotinine	2904	39.4	1.4

Note: (-) Cotinine equivalents represent 50% B/B₀ assay displacement in synthetic matrix.

The compounds having cross-reactivity below 0.39% did not show any significant reaction up to 10 μ g/mL.

ALL THE FOLLOWING HAVE A CROSS-REACTIVITY <0.39%.

Acetaminophen; Acetopromazine; Acetylsalicylic Acid; Amitriptyline; Ascorbic Acid; Benzoic Acid; Caffeine; Chlordiazepoxide; Chlorpromazine; Clenbuterol; Codeine; Dexamethasone; Dextromethorphan; Diclofenac; Dimethyl Sulfoxide; Dipyrone; Doxepin; (-) Ephedrine; Erythromycin; Ethyl-p-amino benzoate; ϵ -amino-n-caproic Acid; Fenoprofen; Flunixin; Folic Acid; Folinic acid; Furosemide; Gemfibrozil; Gentisic Acid; Glipizide; Glutethamide; Glycopyrrolate; Heparin; Hippuric Acid; Hordenine; Hydrocortisone; Ibuprofen; Imipramine; Isoxsuprine; L-Glutamic Acid; Lidocaine; Meperidine; Metaproterenol; Methadone; Methaqualone; Methocarbamol; Methylene Blue; Nalorphine; Naproxen; Niacinamide; Nicotine; Nortriptyline; Orphenadrine; Oxyphenbutazone; PCP; Penicillin G-Potassium; Penicillin G-Procaïne; Pentoxifylline; Phenothiazine; Polyethylene glycol; Prednisolone; Primadone; Procainamide; Promazine; S,S (+) Pseudoephedrine; Pyrantel; Pyrimethamine; Pyrimethamine; Salbutamol; Salicylamide; Salicylic Acid; Sodium Azide; Theophylline; Thiamine; Trimethoprim; Uric Acid

RESULTS INTERPRETATION

Each laboratory should determine the cutoff level for their individual application. When possible, cutoff calibrators and/or standards should be prepared in the same matrix being tested.

Positive Result: Samples with an absorbance less than or equal to the laboratory's designated cutoff calibrator should be considered positive. All positive samples should be confirmed by a quantitative method such as GC/MS.

Negative Result: Samples with an absorbance greater than the laboratory's designated cutoff calibrator should be considered negative.

Qualitative QC Controls: The Neogen positive and negative controls provided in the kit are for QC purposes only. The sole purpose of these controls is to verify that the test kit is performing properly. The controls are not intended for use as cutoff calibrators. The positive control is spiked at 25 ng/mL.

Note: The kit was designed for screening purposes only. It is recommended that all suspect samples be confirmed by a quantitative method such as GC/MS or HPLC.

TECHNICAL SUPPORT

For technical assistance, please contact our technical services department at (859) 254-1221 or email at techservice-toxicology@neogen.com. Representatives are available Monday – Friday from 8:00 am – 6:00 pm EST.

COPYRIGHT

All rights reserved worldwide. No part of this publication may be reproduced, transmitted, transcribed, or stored in any information retrieval system, or translated into any human or computer language in any form or by any means (manual, electronic, mechanical, magnetic, optical, chemical, or otherwise) without expressed written permission.

WARRANTY

Neogen Corporation makes no warranty of any kind, either expressed or implied, except that the material from which its products are made are of standard quality. If any materials are defective, Neogen Corporation will provide a replacement product. Buyer assumes all risk and liability resulting from the use of this product and any of the predictive models. There is no warranty of merchantability of this product for any purpose. Neogen Corporation shall not be liable for any damages, including special or consequential damage, or expense arising directly or indirectly from use of this product.

Products labeled **FORENSIC USE ONLY** are tests that have not been FDA approved.