




Prima[®] Mist

Product Numbers: 353597, 372840, 364876

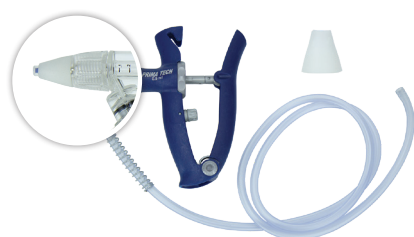

*Prima Mist Sprayer and Replacement Rubber Tip
Syringe sold separately*

Prima Mist	Product Number
Sprayer	353597
Luer Lock Sprayer	372840
Replacement Rubber Tip 12 pk	364876

Key Features

- Intranasal and topical applications
- Molded from chemically resistant clear Tritan[™] material
- Wide dispersion rate with an average particle size of 50-65 microns
- Independent lab tests show a superior product in performance and durability

Product Description

The Prima Mist sprayer is a great addition to our line of syringe accessories. It is engineered for tough multi-use applications. The sprayer screws on to any Prima vaccinator or BMV; it is not suggested for use with syringes over 2 mL or the fast fit system on Prima BMV's. Applications range from intranasal vaccinations to topical antiseptics and anesthetics.

How to Clean Prima Mist Adapter

- Remove adapter from syringe
- Submerge in a solution comprised of dish washing liquid and warm water
- Rinse with water and allow to dry
- It is recommended to replace the rubber tip after cleaning

How to Disinfect Prima Mist Adapter

- As needed spray generously with an aerosolized disinfectant
- Let dry completely

When to Change Prima Mist

- In between pens or in between herds

Syringes that can be used with the Prima Mist

- Line Feed Vaccinators
- Bottle Mount Vaccinators with screw on bottle collar
- Disposable luer lock syringes

Syringes not to use with the Prima Mist

- Bottle Mount Vaccinators with fast fit connectors - this syringe set up has the potential to backfill into the bottle or create enough back pressure to expel the bottle from the syringe

Caution

Do not use the Prima Mist adapter with syringes over 2 mL or with a Fast-Fit connector.

To learn more, call 800.621.8829 (USA/Canada) or 859.254.1221, or visit NEOGEN.com.

