

VERATOX[®] FOR SESAME ALLERGEN



PLEASE READ KIT INSTRUCTIONS COMPLETELY BEFORE PERFORMING TEST.

Product Specifications

Format: Sandwich ELISA

Range of quantification: 2.5-25 ppm sesame

Controls provided: 0, 2.5, 5, 10, and 25 ppm sesame


Testing time: 30 minutes

Product Number: 8530

Call 800.234.5333 to order or visit NEOGEN.com



1 Add 150 μ L controls and extracted samples to red transfer wells.

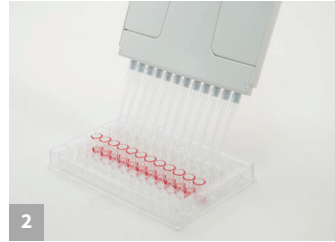

2 Transfer 100 μ L to the antibody wells. Incubate for 10 minutes. Mix for 20 seconds by sliding back and forth on a flat surface.



3 Dump liquid from antibody wells.



4 Wash wells thoroughly. Repeat wash step 5 times with washing solution.



5 Tap out washing solution on absorbent paper towel.



6 Transfer 100 μ L conjugate from reagent boat to antibody wells using a 12-channel pipettor. Incubate for 10 minutes. Mix for 20 seconds by sliding back and forth on a flat surface.



7 Repeat steps 3-5 by dumping out the liquid, thoroughly washing the wells and tapping dry.



8 Transfer 100 μ L substrate from reagent boat to antibody wells using 12-channel pipettor. Incubate for 10 minutes. Mix for 20 seconds by sliding back and forth on a flat surface.



9 Transfer 100 μ L Red Stop from reagent boat to antibody wells. Mix.



10 Read results in a microwell reader.

